


Instinct™


ENDOSCOPIC HEMOCLIP

System Prep 1

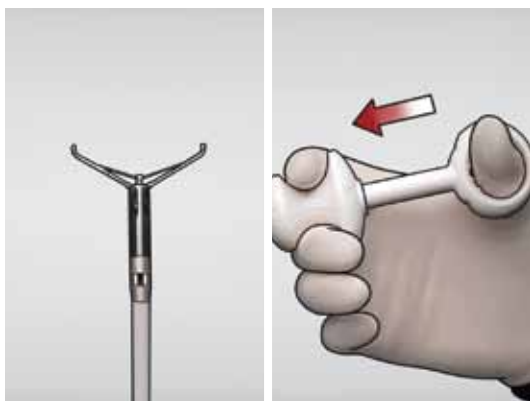

Uncoil device. Verify smooth handle operation and clip action. Open clip by gently moving handle spool distally (away from handle thumb ring). Once clip is fully open, do not continue advancing handle spool as clip may prematurely detach from catheter.


System Prep 2

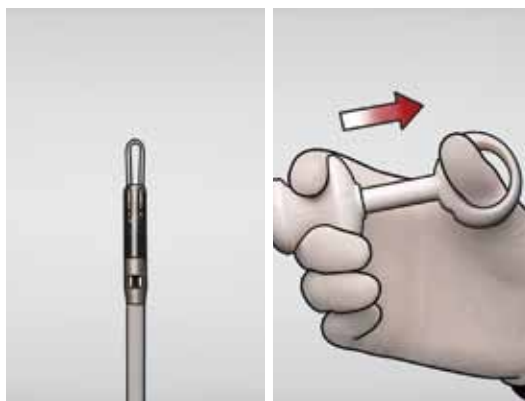

Close clip by moving handle spool proximally until clip is fully closed.

Caution: Do not continue to pull handle spool beyond tactile resistance as this may prematurely deploy clip. Ensure clip is in closed position.

Step 1

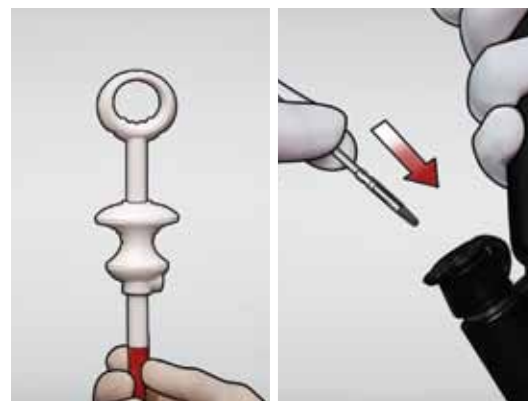

Visually determine desired tissue site. With clip closed and without holding handle spool, advance device in small increments into accessory channel of endoscope.


Caution: Holding handle spool during clip advancement may prematurely deploy clip.

Step 2




Endoscopically confirm that device has exited endoscope. Position distal tip of device toward targeted site.

Step 3

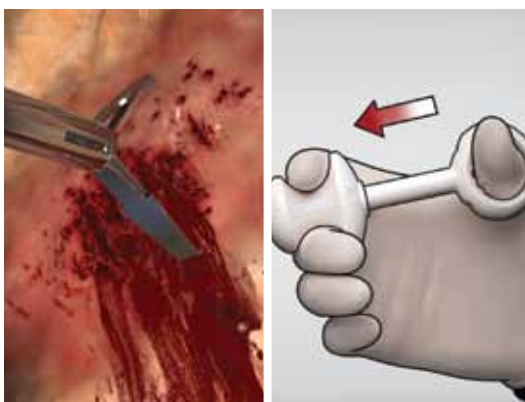

Open clip and advance device into contact with the targeted site.


Step 4

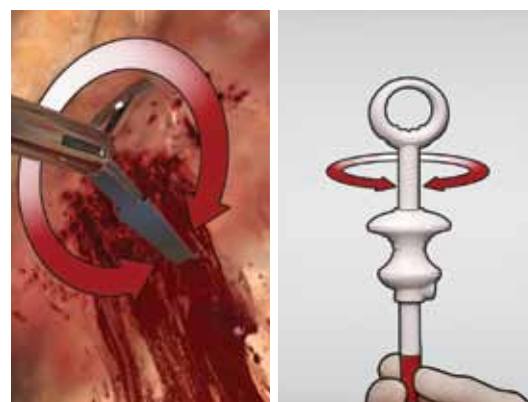

Clip can be rotated by turning handle until clip is in proper position.

Note: Rotation may be limited by clinical circumstances and patient anatomy, among other factors.

Step 5

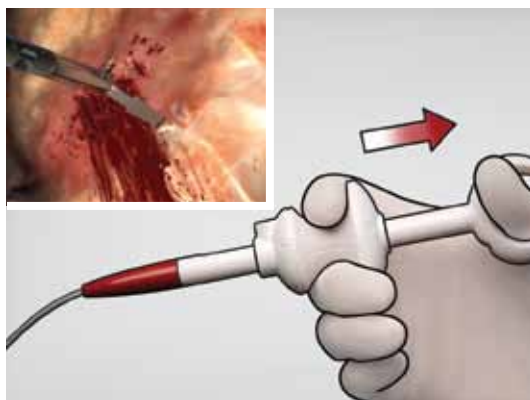


When satisfied with clip position, close clip onto tissue by using slight pressure on handle spool, until tactile resistance is felt.

Note: Ensure clip is pressed firmly against target site for maximum tissue capture.

Caution: Do not continue to pull handle spool beyond tactile resistance until ready to deploy clip; otherwise, clip may not reopen.


Step 6

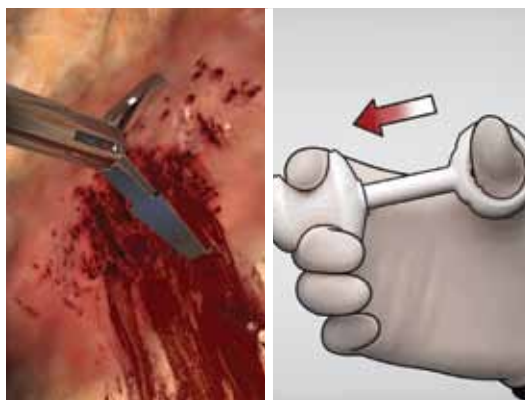

Clip position may now be assessed prior to deployment. If clip is not in desired position, clip may be reopened and repositioned up to 5 times.

Note: Reopening and closing capability may be limited by clinical circumstances and patient anatomy, among other factors.

Step 7

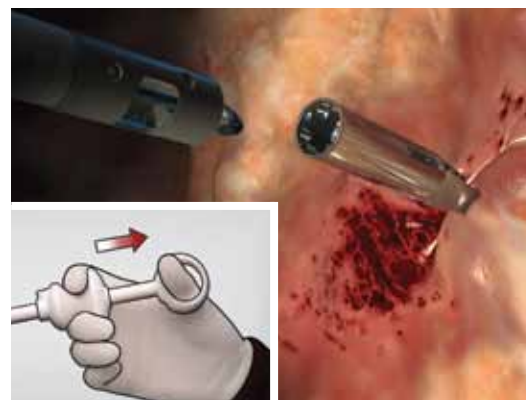


To permanently deploy clip, pull handle spool toward handle thumb ring until clip detaches.

Note: If separation of clip is not immediate, gently move catheter back and forth or use other endoscopic maneuvers to separate catheter from clip.

Caution: Failure to completely release clip may cause clip to pull away from closure site.

Step 8

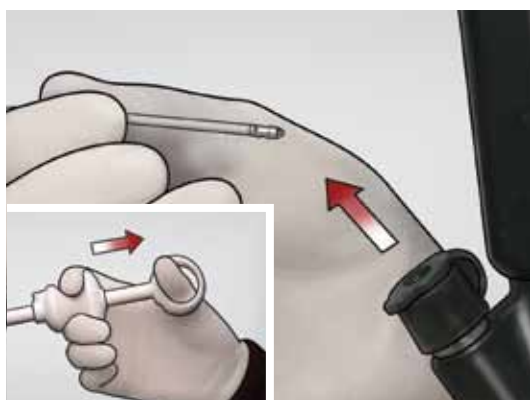

Refer to current instructions for detailed system use.

After clip deployment, continue to apply **slight** pressure on handle spool as device is removed from endoscope.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: orders@cookmedical.com

Australia:
Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283
E-mail: cau.custserv@cookmedical.com

