

Zenith Alpha™

THORACIC ENDOVASCULAR GRAFT

DEVICE PLANNING AND SIZING WORKSHEET


Date: _____ Patient ID: _____
 Hospital: _____
 Physician name: _____
 Physician phone #: _____
 Physician e-mail: _____
 Physician signature: _____
 Purchase order number: _____

Anatomical measurements

Introduction site

- ☐ Right iliac External iliac (EI) measurement _____ mm
☐ Left iliac External iliac (EI) measurement _____ mm

Introduction system: 18-30 mm grafts = 6.0 mm (OD)
 32-38 mm grafts = 7.1 mm (OD)
 40-46 mm grafts = 7.7 mm (OD)

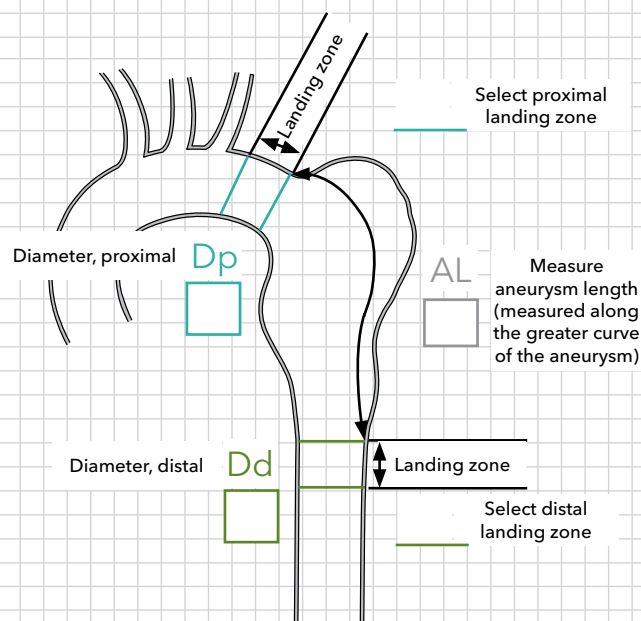

Aortic vessel diameters; CT diameters are outer wall to outer wall

Components

Proximal and distal components are recommended because they provide active fixation proximally and distally. A proximal component should only be used independently for ulcer and/or saccular aneurysms.

Device overlapping

A minimum overlap of 3 stents is recommended. Ensure that the distal graft does not extend beyond the top sealing stent of the proximal graft.

Tapered anatomy

If **Dp** (diameter, proximal neck) is 4 mm > **Dd** (diameter, distal neck), then a proximal tapered graft may be used.

Proximal grafts

Proximal aortic vessel diameter (Dp) =

Select the graft diameter from the vessel oversizing chart below.

plus =
 AL Landing zone Proximal working length

Select the graft length that is closest to the working length.

Proximal Graft (ZTA-P)

Diameter mm	Length mm
18, 20, 22, 24, 26	105
28, 30, 32	109, 155, 201
34, 36	113, 161, 209
38, 40	117, 167, 217
42	121, 173, 225
44, 46	125, 179, 233

Proximal Tapered Graft (ZTA-PT)

Diameter mm	Length mm
26	105
30	108
32	178, 201
34, 36	161, 209
38, 40	167, 217
42	173, 225
44, 46	179, 233

☐ Proximal graft order number = ZTA-P - - -W
 Graft diameter Graft length
☐ Proximal tapered graft order number = ZTA-PT - - - -W
 Graft diameter Graft diameter Graft length

Distal grafts

Distal aortic vessel diameter (Dd) =

Select the graft diameter from the vessel oversizing chart below.

plus =
 AL Landing zone Distal working length

Select the graft length that is closest to the working length.

Distal Graft (ZTA-D)

Diameter mm	Length mm
28, 30, 32	160
34, 36	142, 190
38	147, 197
40	197
42	204
44, 46	211

Distal graft order number = ZTA -D - - -W
 Graft diameter Graft length

Vessel oversizing chart (10%-25%)

Aortic Vessel Outer Diameter (Dp and Dd) mm	Graft Diameter mm
15 - 16	18
17	20
18 - 19	22
20 - 21	24
22 - 23	26
24 - 25	28
26 - 27	30
28 - 29	32
30	34
31 - 32	36
33 - 34	38
35 - 36	40
37 - 38	42
39	44
40 - 42	46

When selecting a graft diameter, refer to the oversizing chart.


Additional considerations may affect the choice of diameter.



www.cookmedical.com

Patient ID: _____

CASE NOTES

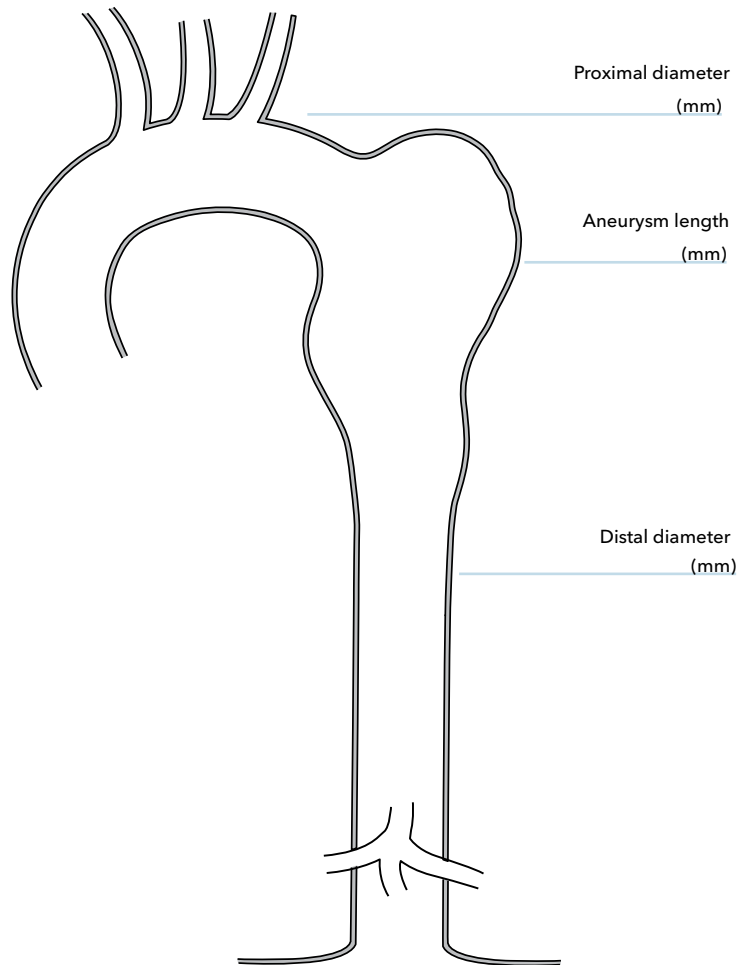

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: orders@cookmedical.com

Australia:
Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
E-mail: cau.custserv@cookmedical.com

