


# Use the Vascular Retrieval Forceps to grasp and retrieve foreign bodies with ease.

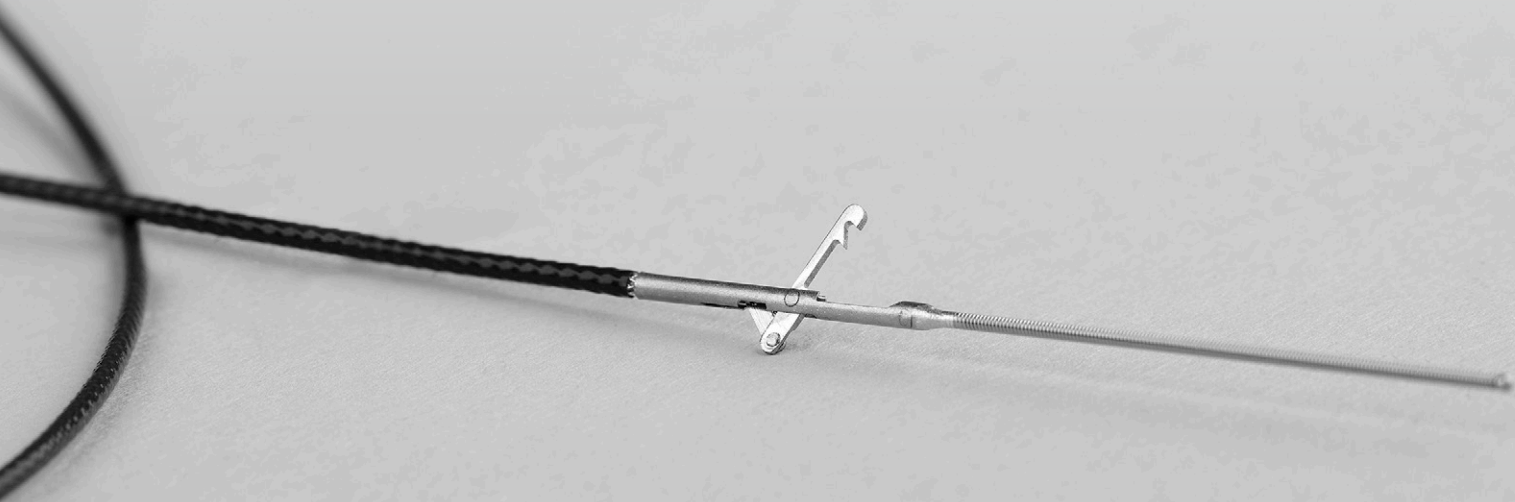

## VASCULAR RETRIEVAL FORCEPS

### **Reach foreign bodies with ease.**

The streamlined design with a straight spring-coiled 2 cm distal tip helps you to secure a firm grasp of the foreign body.

### **Operate simple gripping jaws.**

The easy-to-manipulate jaws lock onto the foreign body, allowing it to be removed from the vessel.

### **View easily under fluoroscopy.**

The braided shaft, jaws, and the distal spring-coil tip are visible under fluoroscopy.





## VASCULAR RETRIEVAL FORCEPS

- The low-profile device can be delivered through a catheter with an ID larger than .045 inches and inserted through a hemostasis valve.
- The device features a PTFE-sheathed, braided stainless steel shaft and strong gripping jaws.
- The distal spring-coil tip is designed to prevent you from inadvertently engaging the vascular wall during manipulations.

Order Number	Reference Part Number	Fr	Length cm
G07187	VRF-3.0-120	3	120

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

### Customer Service

**EMEA:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)

Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)

Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)

Belgium: +32 27001633, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)

Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)

Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)

France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)

Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)

Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)

Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)

Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)

Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)

Norway: +47 23162968, [no.orders@cookmedical.com](mailto:no.orders@cookmedical.com)

Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)

Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)

Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)

Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)

Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)

United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)

[www.cookmedical.com](http://www.cookmedical.com)

**Americas:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: [orders@cookmedical.com](mailto:orders@cookmedical.com)

**Australia:**

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)

