


# Cross lesions with just one set.

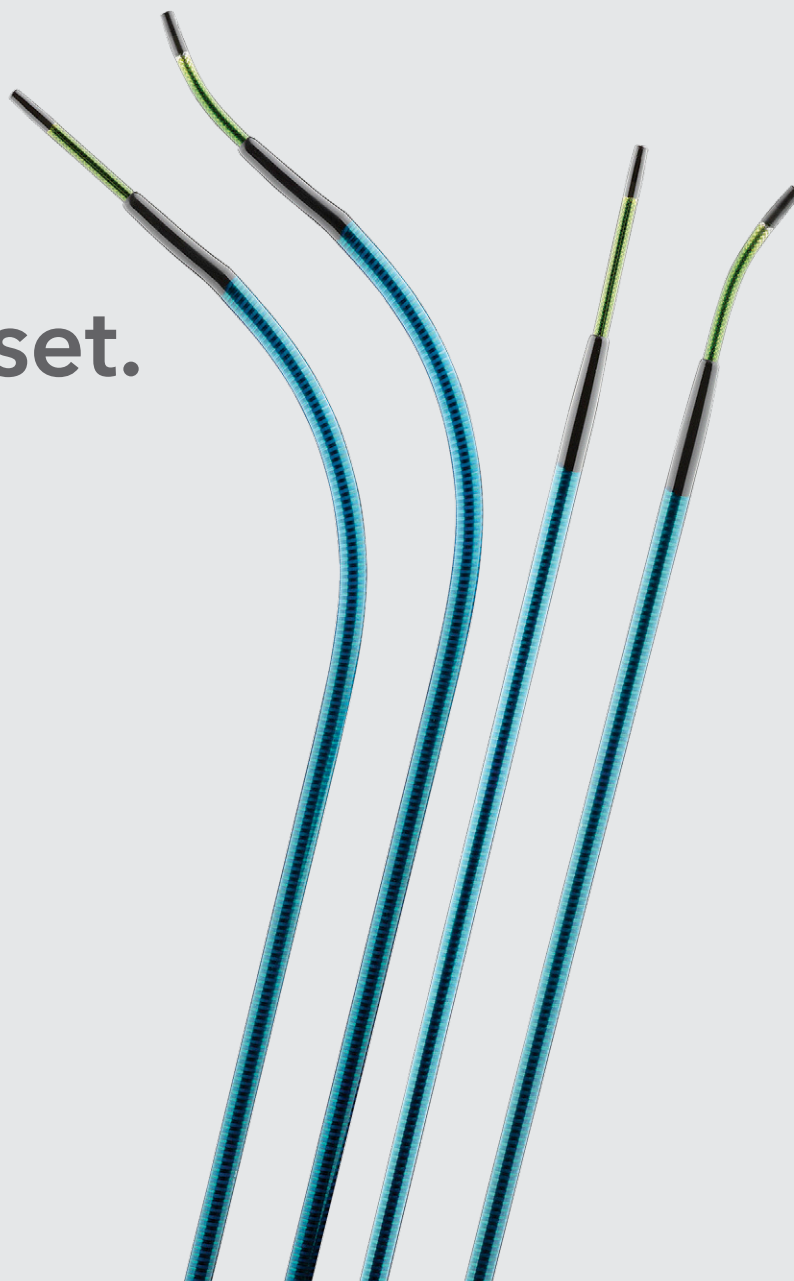

## TriForce® PERIPHERAL CROSSING SET

### Enhanced penetration

The enhanced crossing ability of the CXI® catheter joins with the pushability of the Flexor® sheath to cross vascular lesions.

### Strong pushability

A tungsten-loaded tip aids pushability and adds to radiopacity as it goes through the lesion. The hydrophilic coating on both the CXI and the Flexor sheath enhances trackability.


### Tapered transition

The 5 Fr sheath tapers to a 4 Fr catheter that hugs a .035 inch diameter wire guide for an uninterrupted transition.

Data on file



[www.cookmedical.com](http://www.cookmedical.com)

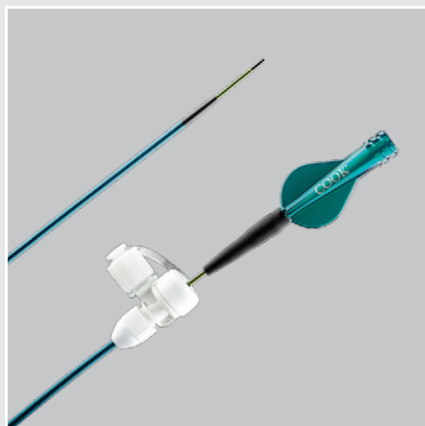

## TriForce®

### PERIPHERAL CROSSING SET

- TriForce Peripheral Crossing Set includes an outer 5.0 Fr Flexor® sheath, a 4.0 Fr CXI® Support Catheter, and a Peel-Away® sheath.
- Flexor outer sheath tapers to a 4.0 Fr inner diameter at the tip to ensure a smooth transition to the CXI Support Catheter.
- Flexor outer sheath features radiopaque tip and hydrophilic coating.
- Accepts a .035 inch diameter wire guide.
- The infusion dynamic pressure of the CXI Support Catheter should not exceed 1,200 psi (82.7 bar). The static pressure should not exceed 700 psi (48.3 bar).

Order Number	Reference Part Number	Catheter Fr/Length cm	Sheath ID/Length Fr/cm	Sheath Tip Fr	Accepts Wire Guide Diameter Inch
<b>Set with Straight Tip Sheath and Straight Tip Catheter</b>					
G56412	KCXS-5.0-35-65-RB-0/0-HC	4.0/65	5.0/55	4.0	.035
G56416	KCXS-5.0-35-100-RB-0/0-HC	4.0/100	5.0/90	4.0	.035
<b>Set with Straight Tip Sheath and Curved Tip Catheter</b>					
G56413	KCXS-5.0-35-65-RB-0/DAV-HC	4.0/65	5.0/55	4.0	.035
G56417	KCXS-5.0-35-100-RB-0/DAV-HC	4.0/100	5.0/90	4.0	.035
<b>Set with Curved Tip Sheath and Straight Tip Catheter</b>					
G56414	KCXS-5.0-35-65-RB-MPB/0-HC	4.0/65	5.0/55	4.0	.035
G56418	KCXS-5.0-35-100-RB-MPB/0-HC	4.0/100	5.0/90	4.0	.035
<b>Set with Curved Tip Sheath and Curved Tip Catheter</b>					
G56415	KCXS-5.0-35-65-RB-MPB/DAV-HC	4.0/65	5.0/55	4.0	.035
G56419	KCXS-5.0-35-100-RB-MPB/DAV-HC	4.0/100	5.0/90	4.0	.035

Set includes 5.0 Fr Flexor outer sheath, 4.0 Fr CXI Support Catheter and Peel-Away sheath.

Not all part numbers shown on this product information sheet may be approved for sale in all regulatory jurisdictions. Consult with your local Cook representative or customer service center for details.

### Customer Service

**EMEA:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)  
Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)  
Belgium: +32 27001633, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)  
Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)  
Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)  
France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)  
Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)  
Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)  
Norway: +47 23162968, [no.orders@cookmedical.com](mailto:no.orders@cookmedical.com)  
Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)  
Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)  
Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)

[www.cookmedical.com](http://www.cookmedical.com)

**Americas:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335  
E-mail: [orders@cookmedical.com](mailto:orders@cookmedical.com)

**Australia:**  
Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283  
E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)

