


Get glowing results.

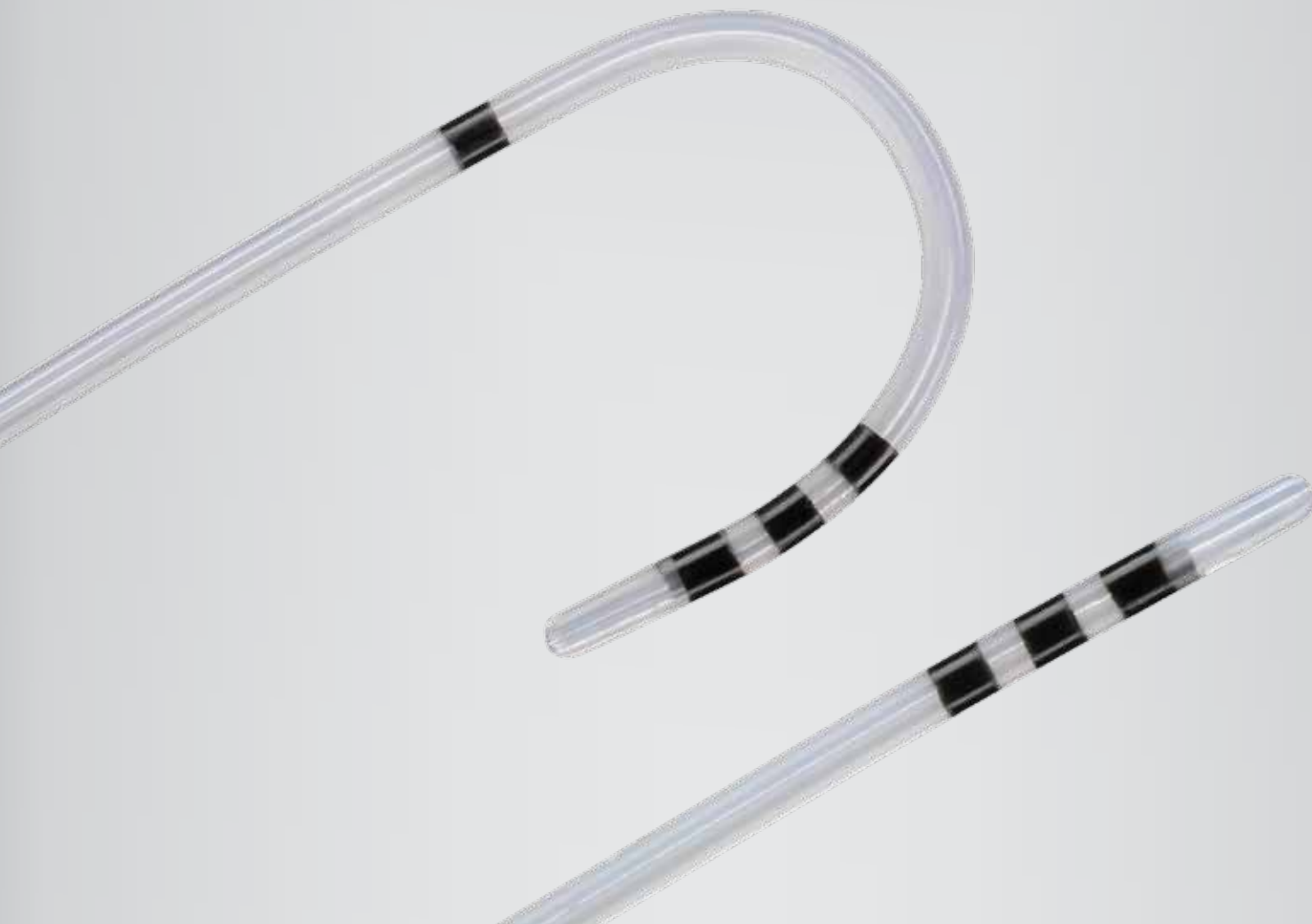

Glo-Tip II®

DOUBLE LUMEN ERCP CATHETER

The Cook Glo-Tip II catheter balances two significant factors in cannulating procedures – access and visualization. The catheter's DomeTip® permits smooth cannulation with potentially less trauma and a radiopaque band inside the wire lumen emits a distinctive "glow" under fluoroscopy. And, for convenience, its double lumen design permits contrast injection through a dedicated lumen while the catheter is placed over a wire guide. The Glo-Tip II: part of Cook's full line of innovative endoscopic accessories designed to meet your patients' needs.



www.cookmedical.com



Glo-Tip II®

DOUBLE LUMEN ERCP CATHETER

This device is used for endoscopic cannulation of the ductal system.

Supplied sterile and is disposable – intended for single use only.

- DomeTip® provides smoother access with potentially less trauma
- Radiopaque tip provides excellent fluoroscopic visualization
- Double lumen design permits contrast injection through a dedicated lumen while catheter is placed over a wire guide
- Mark V System® – with endoscopic markers at 5, 7, 9 and 12 cm – can be used to determine stent length and to verify cannulation depth

Order Number	Reference Part Number	Catheter Size Fr	Catheter Length cm	Catheter Tip cm	Tip	Remarks	Color Connection
G22226	GT-2-T	6	200	100	Dometip	Standard	Purple
G24694	GT-2-T-RB	6	200	100	Dometip	Radiopaque Bands	Purple
G24506	GT-2-T-ANG	6	200	100	Dometip	Angled	Purple

Minimum accessory channel 2.8 mm. Stylet length 100 cm. Recommend Tracer Metro® Direct™ .035 inch wire guide, sold separately.

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do

Distributors: +353 61239240, ssc.distributors@cookmedical.com

Austria: +43 179567121, oe.orders@cookmedical.com

Belgium: +32 27001633, be.orders@cookmedical.com

Denmark: +45 38487607, da.orders@cookmedical.com

Finland: +358 97251996, fi.orders@cookmedical.com

France: +33 171230269, fr.orders@cookmedical.com

Germany: +49 6950072804, de.orders@cookmedical.com

Hungary: +36 17779199, hu.orders@cookmedical.com

Ireland: +353 61239252, ie.orders@cookmedical.com

Italy: +39 0269682853, it.orders@cookmedical.com

Netherlands: +31 202013367, nl.orders@cookmedical.com

Norway: +47 23162968, no.orders@cookmedical.com

Spain: +34 912702691, es.orders@cookmedical.com

Sweden: +46 858769468, se.orders@cookmedical.com

Switzerland - French: +41 448009609, fr.orders@cookmedical.com

Switzerland - Italian: +41 448009609, it.orders@cookmedical.com

Switzerland - German: +41 448009609, de.orders@cookmedical.com

United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283

E-mail: cau.custserv@cookmedical.com

