


Increase your procedural accuracy
with balloon-expandable stent
technology engineered for precision.



Formula 418[®] Biliary

BALLOON-EXPANDABLE STENT

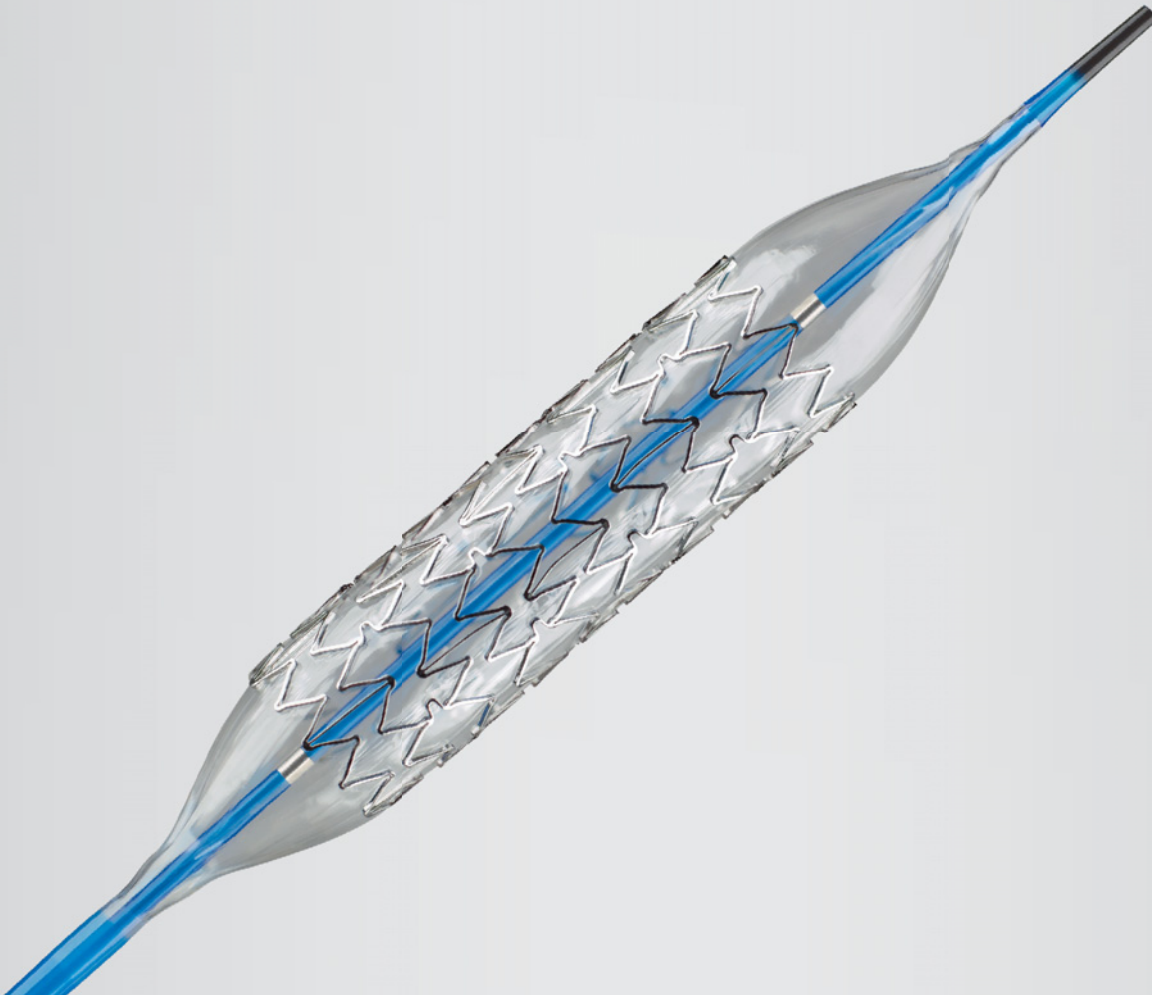

www.cookmedical.com

When precise deployment matters, so does technology.

To precisely place most balloon-expandable stents, you have to be able to predict how much they will shorten. But precision shouldn't involve guesswork. That's why Formula stents are engineered for precision: precise placement, precise treatment and precise conformability.

Precise Placement


Formula stents don't shorten, and their ends align with the balloon markers. These features help you to place the stent and to treat the stricture reliably and precisely.

Precise Conformability


Formula balloon-expandable stents have additional crowns that give their stents a more circular shape and improve their scaffolding.

Precise Treatment

With one of the lowest tip profiles on the market, Formula 418 Biliary can help you cross strictures. The minimal balloon overhang facilitates precise treatment.

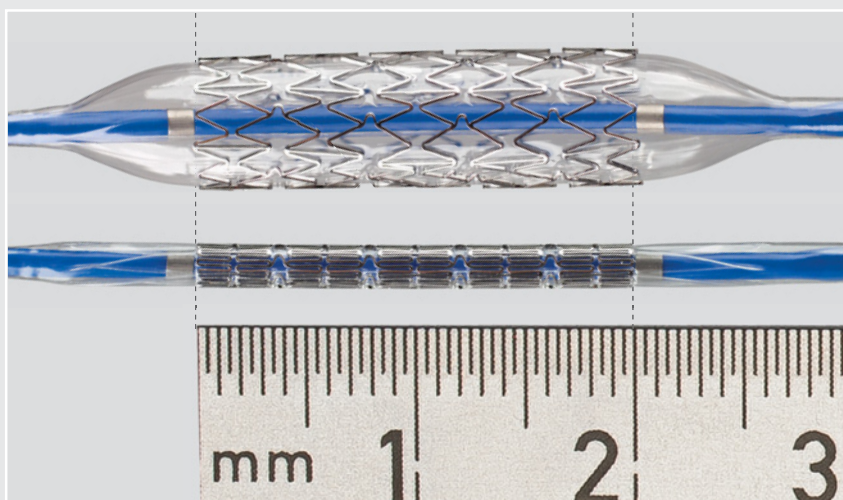


Precise Placement

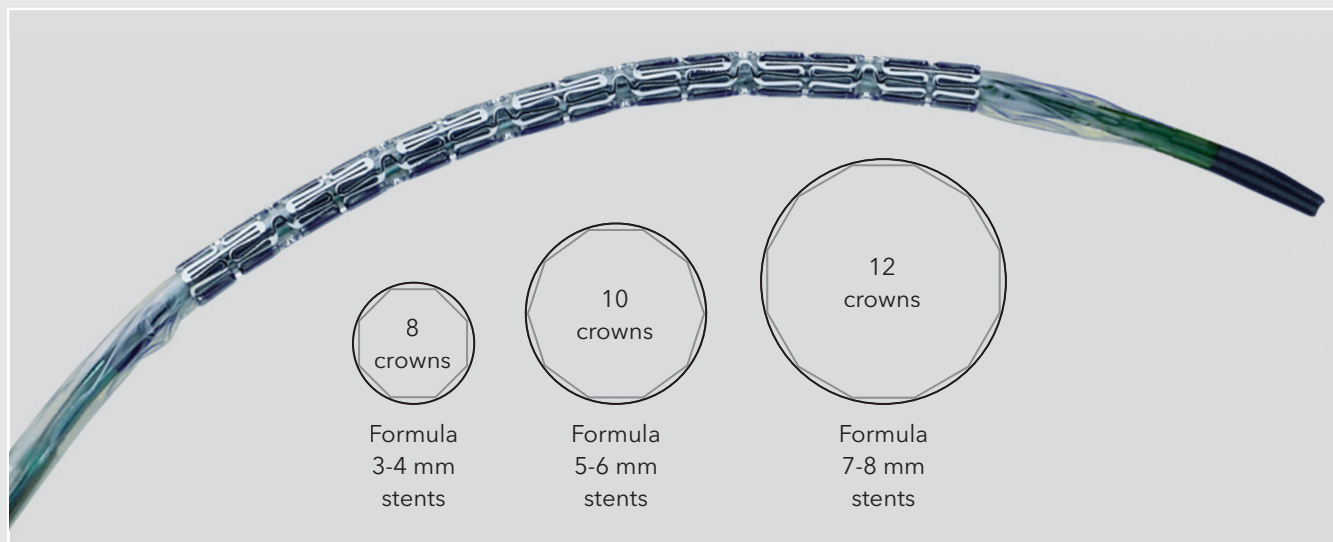

Precise Conformability

Order Number	Reference Part Number	Stent Diameter mm	Stent Length mm	Delivery System Length cm	Sheath/ Guiding Catheter Fr	Nominal Pressure atm	Rated Burst Pressure atm
Formula 418 – 80 cm Over-the-Wire Delivery System							
G42914	FORB4-18-80-3-12	3	12	80	5.0/6.0	8	12
G42915	FORB4-18-80-3-16	3	16	80	5.0/6.0	8	12
G42916	FORB4-18-80-3-20	3	20	80	5.0/6.0	8	12
G42917	FORB4-18-80-3-24	3	24	80	5.0/6.0	8	12
G42918	FORB4-18-80-3-30	3	30	80	5.0/6.0	8	12
G42922	FORB4-18-80-4-24	4	24	80	5.0/6.0	8	12
G42923	FORB4-18-80-4-30	4	30	80	5.0/6.0	8	12
G42927	FORB4-18-80-5-24	5	24	80	5.0/6.0	8	12
G42928	FORB4-18-80-5-30	5	30	80	5.0/6.0	8	12
G42932	FORB4-18-80-6-24	6	24	80	5.0/6.0	8	12
G42933	FORB4-18-80-6-30	6	30	80	5.0/6.0	8	12
G42937	FORB4-18-80-7-24	7	24	80	6.0/7.0	8	12
G42938	FORB4-18-80-7-30	7	30	80	6.0/7.0	8	12
G42939	FORB4-18-80-8-12	8	12	80	6.0/8.0	8	12
G42940	FORB4-18-80-8-16	8	16	80	6.0/8.0	8	12
G42941	FORB4-18-80-8-20	8	20	80	6.0/8.0	8	12
G42942	FORB4-18-80-8-24	8	24	80	6.0/8.0	8	12
G42943	FORB4-18-80-8-30	8	30	80	6.0/8.0	8	12

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

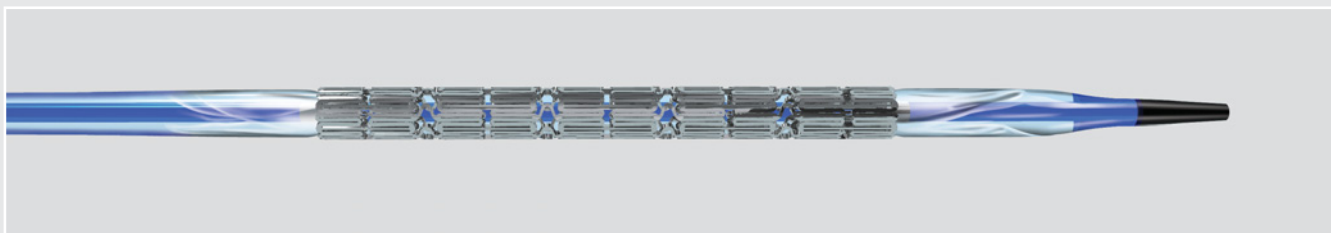

Precise Treatment

Order Number	Reference Part Number	Stent Diameter mm	Stent Length mm	Delivery System Length cm	Sheath/ Guiding Catheter Fr	Nominal Pressure atm	Rated Burst Pressure atm
Formula 418 – 135 cm Over-the-Wire Delivery System							
G42944	FORB4-18-135-3-12	3	12	135	5.0/6.0	8	12
G42945	FORB4-18-135-3-16	3	16	135	5.0/6.0	8	12
G42946	FORB4-18-135-3-20	3	20	135	5.0/6.0	8	12
G42947	FORB4-18-135-3-24	3	24	135	5.0/6.0	8	12
G42948	FORB4-18-135-3-30	3	30	135	5.0/6.0	8	12
G42952	FORB4-18-135-4-24	4	24	135	5.0/6.0	8	12
G42953	FORB4-18-135-4-30	4	30	135	5.0/6.0	8	12
G42957	FORB4-18-135-5-24	5	24	135	5.0/6.0	8	12
G42958	FORB4-18-135-5-30	5	30	135	5.0/6.0	8	12
G42962	FORB4-18-135-6-24	6	24	135	5.0/6.0	8	12
G42963	FORB4-18-135-6-30	6	30	135	5.0/6.0	8	12
G42967	FORB4-18-135-7-24	7	24	135	6.0/7.0	8	12
G42968	FORB4-18-135-7-30	7	30	135	6.0/7.0	8	12
G42969	FORB4-18-135-8-12	8	12	135	6.0/8.0	8	12
G42970	FORB4-18-135-8-16	8	16	135	6.0/8.0	8	12
G42971	FORB4-18-135-8-20	8	20	135	6.0/8.0	8	12
G42972	FORB4-18-135-8-24	8	24	135	6.0/8.0	8	12
G42973	FORB4-18-135-8-30	8	30	135	6.0/8.0	8	12

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: orders@cookmedical.com

Australia:
Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283
E-mail: cau.custserv@cookmedical.com

