


Light your way.



Bush

URETERAL ILLUMINATING CATHETER SET

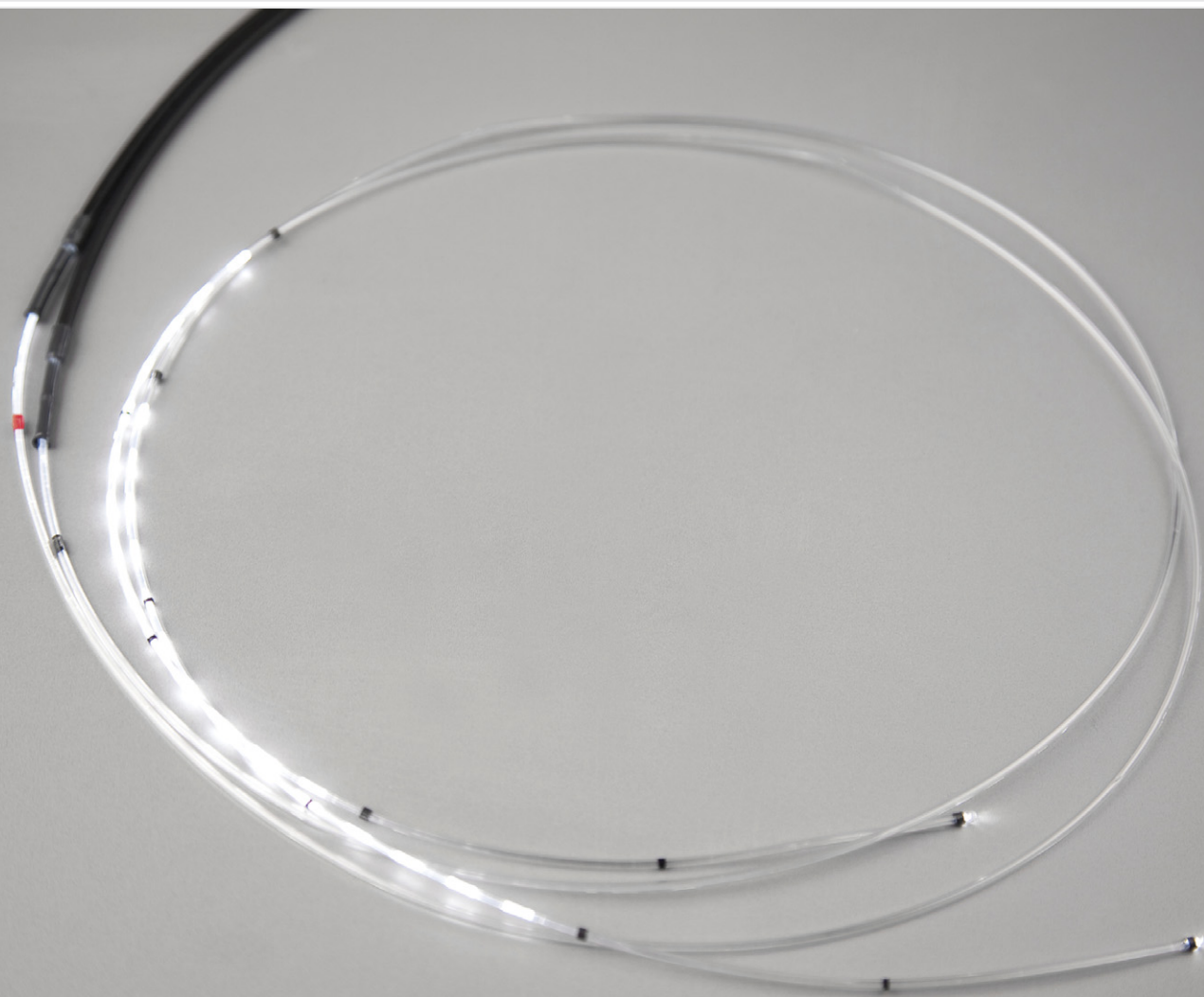

www.cookmedical.com

Bush

URETERAL ILLUMINATING CATHETER SET

When connected to a light source, the Bush Ureteral Illuminating Catheter Set allows you to visualize the ureters during surgery.

The Bush catheters can be used with a variety of light sources, so a dedicated light source is not required. Choose between two catheter configurations to make bilateral stenting possible while using either a shared light source or a separate light source for each ureter.

The Bush catheters are available in either a single- or dual-lumen set. The dual-lumen catheter is designed for drainage during surgery.

Features

Light fiber extensions

A special fitting allows you to disconnect the catheter to allow for scope removal after catheter placement.

Color-coded proximal ends

The color coding allows for easy external identification of the left and right ureters.

Radiopaque marks

The radiopaque marks ensure proper placement of the catheter within the ureter.




Bush SL illuminating a ureter during laparoscopic hysterectomy.


Notes on use:

- Keep the light source on when in use. Dim or turn off the light when not in use.
- Always turn on the light source at the minimum setting and increase the light intensity as necessary.
- Dim the room and reduce laparoscopic lighting for better visualization.
- Single or dual light sources may be needed, depending on your preference and the catheter's configuration. Check the product label to verify whether the Bush catheter is the single-lumen or dual-lumen configuration.
- Contact your sales representative for in-servicing of the device.


The Bush catheters work seamlessly with the VueLite® LED Light Source to produce bright light. The VueLite is a practical alternative to traditional light sources because it uses a low amount of energy and generates a low amount of heat.¹ The portability of VueLite also allows you to perform a procedure without being tethered to a power source.

The catheters are also compatible with many other commonly available light sources. Depending on the light source, an adapter may be required to successfully use the Bush DL or Bush SL catheter.

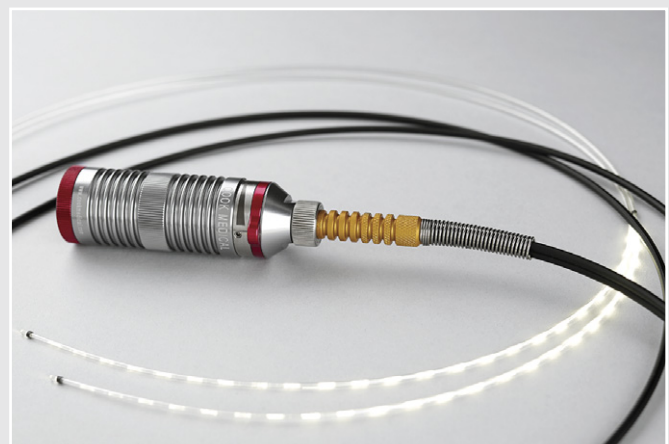

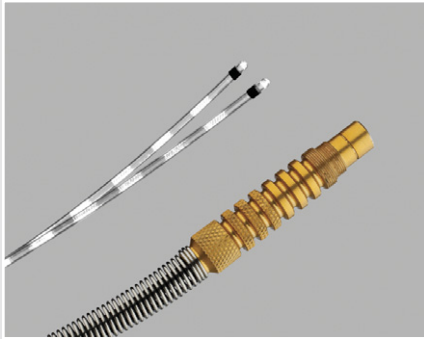

The Bush Ureteral Illuminating Catheter Set is used for transillumination of the ureters during laparoscopic or open surgical procedures. Transillumination helps to identify, and minimize potential for trauma to, the ureters.

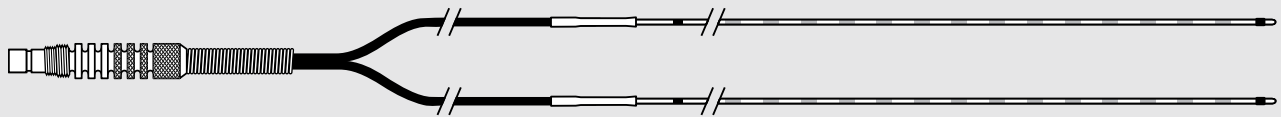

Order Number	Reference Part Number	Fr	Length cm	Shielded Light Fiber Length cm
Two Single-Lumen Catheters				
With Two Light Source Plugs				
G16745	084100	5.0	80	130
With Common Light Source Plug				
G16746	084120	5.0	80	130
Two Dual-Lumen Catheters				
With Two Light Source Plugs				
G16143	084510	7.0	80	130
With Common Light Source Plug				
G16747	084520	7.0	80	130

Note: The light guide plug that is supplied fits ACMI light sources MV9082, MV9083, and ALV-1.

For Olympus light sources, obtain ACMI light source adapter 8914 from ACMI.

For STORZ light sources, obtain ACMI plug adapter 487 A from KARL STORZ.

For STORZ xenon light sources, obtain turret adapter 487UO from KARL STORZ.

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.



Customer Service - Urology

EMEA: EDI - www.cookmedical.com/edi.do
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Poland: +48 223060159, pl.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: orders@cookmedical.com

Australia:
Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
E-mail: cau.custserv@cookmedical.com

